
Specification and Design Guide

Vintage Series

Surface Mounted Vertical Headwall

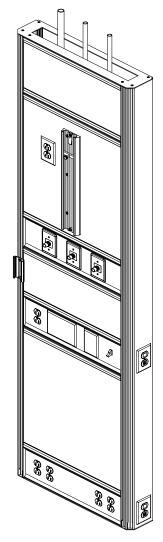

* Unit shown: Vintage Series 22" Full length Headwall

Item number definition

VHW-XXYY "XX"= Width in inches "YY"= Height in inches

Available Height:

Half (05) Three-Quarter (07) Full (10)

Construction Specification

The Headwall Unit shall be constructed of .094" thick 6063-t 5 alloy, extruded aluminium for the structural frame members. All exposed aluminium surface shall be clear anodized. Each side corner shall be equipped with Medical Design's 45° two track system as an integral part of the main frame, designed to manage hospital equipment and accessories both horizontally and vertically. An all aluminium fascia strip shall receive medical gas outlets and electrical devices. The decorative laminated panel shall be constructed of high density fire retardant particle board. This product is cULus listed.

Medical Gas

Medical gas outlets shall be supplied, installed and manifolded by Medical Design. Connection points for copper tubing are located 6" above top of unit. Medical gas outlet will be manifolded with copper piping unless otherwise specified. Installation and testing shall comply with National Fire Protection Association (NFPA 99c) & C.R.N. (BNQ 5710-500, CSA 305.1 / B51-95).

Electrical Equipment

Electrical devices shall be supplied and pre-wired by Medical Design. Wiring shall run through metal raceway and terminate at the top of the unit connected to the Hospital electrical system. Conductors shall be No 12 stranded RW-90 copper wire for normal and emergency power circuits. All devices shall be grounded with #10 AWG green jacket stranded copper wire. All electrical receptacles are NEMA 5-15/20R, 125 volts AC, hospital grade and comply with CSA and are cULus listed.

Emergency receptacle shall be red in color and ivory for normal power.

Communication

Communication devices shall be supplied and installed by others. Provisions for devices shall be coordinated with Medical Design.

Site Installation

All units shall be installed in accordance with Medical Design's installation instructions. The mechanical and electrical contractor shall provide the final connection to the hospital System. The contractor has the responsibility to insure that every unit is subject to a final site inspection in compliance with all National, State and Local Codes.






Typical Shop Drawing



VINTAGE SERIES GENERAL CARE PRIVATE ROOM HEADWALL

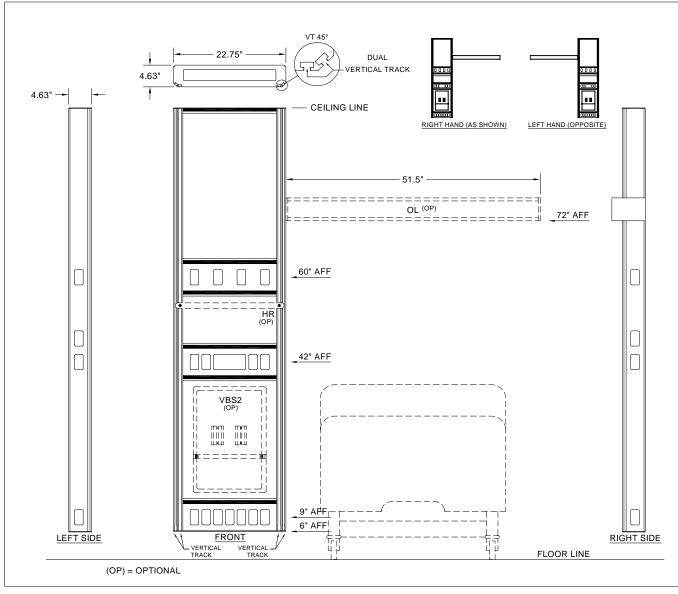

MEDICAL GAS		
ID	QTY	DESCRIPTION
02	-	OXYGEN OUTLET; Q.C. OU D.I.S.S.
Α	-	MEDICAL AIR OUTLET; Q.C. OU D.I.S.S.
V	-	VACUUM OUTLET; Q.C. OU D.I.S.S.
N_2	-	NITROGEN OUTLET; D.I.S.S. ONLY
N_2O	-	NITROUS OXIDE OUTLET; Q.C. OU D.I.S.S.
EVAC	-	EVACUATION OUTLET; Q.C. OU D.I.S.S.
CO ₂	-	CARBON DIOXIDE OUTLET; D.I.S.S. ONLY
	1	ELECTRICITY AND COMMUNICATION
DR	-	DUPLEX RECEPTACLE, 20A, 125V, H.G., (RED) EMERGENCY POWER (HUBBELL #8300HR)
DW	-	DUPLEX RECEPTACLE, 20A, 125V, H.G., (WHITE) NORMAL POWER (HUBBELL #8300HW)
SR	-	SIMPLEX RECEPTACLE, 20A, 125V, H.G., (RED) EMERGENCY POWER (HUBBELL #8310R)
SW	-	SIMPLEX RECEPTACLE, 20A, 125V, H.G., (WHITE) NORMAL POWER (HUBBELL #8310W)
LS	-	LIGHT SWITCH, 1P, 15A, 125V, (WHITE)
TW	-	TWIST-LOCK RECEPTACLE, 20A, 125V, H.G.
GR	-	GROUNDING RECEPTACLE, HAMPDEN SLR-3S
МО	-	MONITOR OUTLET PROVISION, 2 GANGS
NC	-	NURSE CALL PROVISION, 3 GANGS
СВ	-	CODE BLUE PROVISION, 1 GANG
МО	-	MONITOR OUTLET PROVISION, 2 GANGS
D	-	DATA PROVISION, 1 GANG
Т	-	TELEPHONE PROVISION, 1GANG
В	-	BLANK FACEPLATE, 1 GANG, STAINLESS STEEL
ACCESSORIES		
HR	-	HORIZONTAL RAIL
МВ	-	MONITOR BRACKET, 19"
OL	-	OVERBED LIGHT, 4'
SL	-	VACUUM BOTTLE SLIDE
VBS2	-	VACUUM BOTTLE STORAGE COVE (2 BOTTLES) WITH 2 SLIDES AND 1 GUARD RAIL
VT	4	VERTICAL TRACK
OTHER SERVICES AND ACCESSORIES ARE AVAILABLE. CONTACT US FOR MORE INFORMATION.		